
Asguard®

ISD001, ISD002, ISD003, ISD004, ISD006, ISD008 & ISD030

Non-Woven Island Dressing *(Instructions For Use)*


Asguard Non-Woven Island Dressings feature an adhesive border that securely seals the dressing in place, helping to prevent contamination and ensuring the wound remains protected from external factors. This adhesive is skin-friendly and provides reliable fixation, yet is easy to remove, minimising trauma to the healing tissue.

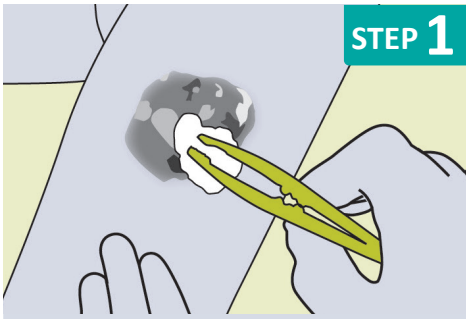


Cleanse wound thoroughly, using local protocols. The skin surrounding the wound should be clean and dry.




Dressing size is dependent on the wound surface area – the dressing pad should extend 1-2cms beyond the wound margin.

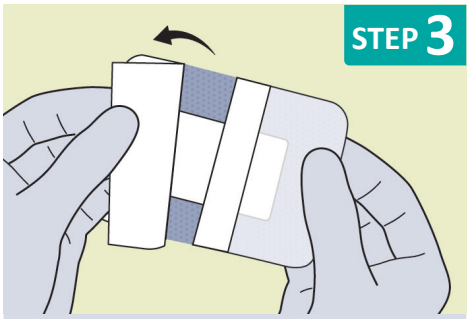


Open the sterile packaging and remove the dressing. Separate the paper backing from the Asguard Non-Woven Island Dressing.

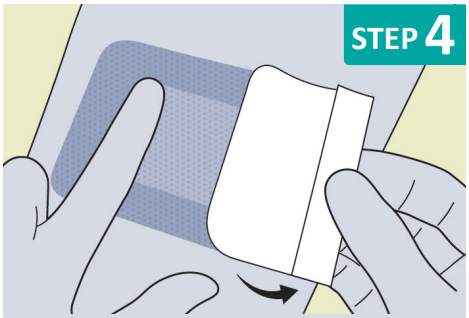

Apply the Asguard Non-Woven Island Dressing with central pad covering entire wound bed. Apply gentle pressure around the edge of pad.

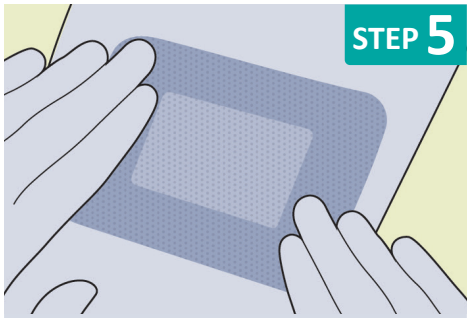

Smooth down the edges of the non-woven border ensuring adherence.

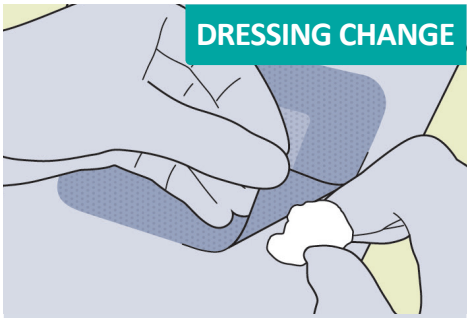

The clinician-in-charge is responsible for determining the need for dressing changes, dependent on the stage and phase of wound healing and exudate level.

The dressing should be changed when exudate absorption is considered to be extending beyond the wound surface area.

7 days is the maximum period between dressing changes.

To minimise potential damage to fragile skin, Asguard Non-Woven Island Dressing should be eased off using an adhesive remover wipe.

INDICATIONS

Asguard Non-Woven Island Dressings are indicated for use as primary dressings in the management of all types of superficial wounds with light exudate and where water resistance is not required, such as:

- Abrasions and lacerations
- Blisters and minor skin loss
- Skin tears
- Post-operative & post-procedural sites

CONTRAINDICATIONS

Should not be used if allergies to any of its components, is known.

GENERAL INFORMATION

If infection is present, a primary wound cleansing dressing e.g. Osmonate Calcium Alginate dressings or rope, or an antimicrobial layer such as Zorflex should be used in combination with Asguard Silicone Hydroporous Foam. The clinician in-charge is responsible for determining if necessary.

At Sentry, we are dedicated to upholding sustainable principles to safeguard our planet, people, and future, and we strive to consistently act responsibly in doing so.

POUCH - PAPER	POUCH - FILM	BOX
RECYCLE ME	GENERAL WASTE	RECYCLE ME

Sentry Medical Pty Ltd
1300 995 999 | sales@sentrymedical.com.au
www.sentrymedical.com.au

Sentry MEDICAL

SMHO-MKT-IFU-0008-REV 2_10_26 Seek advice from a health professional before use. All product information is correct at date of publication. The image may differ slightly from the actual product.