



SureSafe Poly-Reinforced Surgical Gown

SureSafe® poly-reinforced gowns are made from 4-layer polypropylene materials with anti-static and fluid resistance. Two layers of spunbond polypropylene ensure an extraordinary flexibility and drapeability, provide tough abrasion resistance and cloth like feel, and two layers of meltblown polypropylene forming an exceptional microbial barrier against bacteria. The critical zones including arm and chest areas are protected by poly-reinforcement with an additional strong and breathable membrane. The seams on the sleeves are ultrasonic sealed to avoid fluid penetration through pin holes.

Compliance AAM Level 2 - EN 13795 – Full Coverage

4-layer anti-static breathable fabric

Fluid repellent

Excellent microbial barrier

Light weight material, less recycled waste

Velcro neck binding

Generous cut

Large transfer wrap – 90 x 90cm

Cost effective and environmentally friendly

2 paper towels standard or 1 cotton huck towel*

Sterile individually packed, available non-sterile*


*Contact Aaxis sales consultant for lead time and MOQ



PRODUCT CODE	MATERIAL	STERILE REINFORCED		SIZE	COTTON / PAPER TOWEL	AAMI LEVEL Full Coverage	SKU/TU (MOQ)	SKU/ CTN
S61702	35gsm 4 Layer SMMS	No	Yes	Medium	-	Level 2	54	54
S61703	35gsm 4 Layer SMMS	No	Yes	Large	-	Level 2	54	54
S71753	35gsm 4 Layer SMMS	Yes	Yes	Large	1 cotton	Level 2	24	24
S71754	35gsm 4 Layer SMMS	Yes	Yes	X-Large	1 cotton	Level 2	24	24
S71755	35gsm 4 Layer SMMS	Yes	Yes	XX-Large	1 cotton	Level 2	24	24
S71752P	35gsm 4 Layer SMMS	Yes	Yes	Medium	2 paper	Level 2	30	30
S71753P	35gsm 4 Layer SMMS	Yes	Yes	Large	2 paper	Level 2	30	30
S71754P	35gsm 4 Layer SMMS	Yes	Yes	X-Large	2 paper	Level 2	30	30

SureSafe Poly-Reinforced Surgical Gown

Product Dimensions

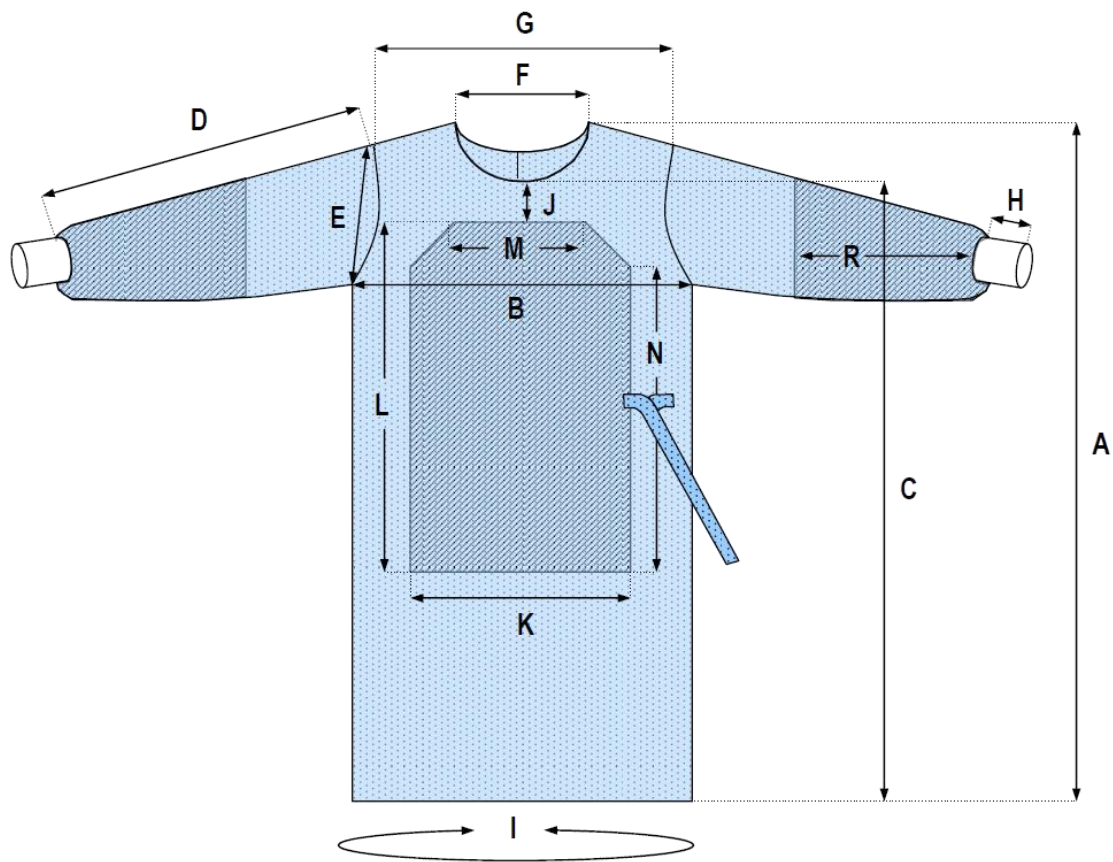

PRODUCT CODE	SIZE	DIMENSIONS														
		A	B	C	D	E	F	G	H	I	J	K	L	M	N	O
S71752*	Medium	120	72	112	55	26	22	57	8	150	9	50	60	32	32	30
S71753*	Large	128	72	118	58	27	24	65	8	156	9	50	90	38	36	-
S71754*	X Large	143	72	134	60	28	24	64	8	160	10	50	90	38	36	-
S71755*	XX Large	148	85	140	65	32	26	77	8	175	12	50	90	42	40	-
	Tolerance	±2	±1.5	±2	±1.5	±1	±1	±1.5	±0.75	±2	±0.75	±1.5	±1.5	±1	±2	

*Sizing applies to any Suresafe 35 Poly-reinforced gown sterile or non sterile, with 1 cotton huck towel or 2 paper towel configurations