



# SureSafe Poly-Reinforced Surgical Gown

SureSafe® poly-reinforced gowns are made from 5-layer polypropylene materials with anti-static and fluid resistance. Two layers of spunbond polypropylene ensure an extraordinary flexibility and drapeability, provide tough abrasion resistance and cloth like feel, and two layers of meltblown polypropylene forming an exceptional microbial barrier against bacteria. The critical zones including arm and chest areas are protected by poly-reinforcement with an additional strong and breathable membrane. The seams on the sleeves are ultrasonic sealed to avoid fluid penetration through pin holes.

5-layer anti-static breathable fabric

Compliance AAM Level 3 - EN 13795 – Full Coverage

Fluid repellent

Excellent microbial barrier

Excellent physical Barrier to contaminants/fluids

Stronger fabric, matte finish, resistance to tears

Light weight material, less recycled waste

Softer fabric with comfort inside and outside

Generous cut

Velcro neck binding

Large transfer wrap – 90 x 90cm

Cost effective and environmentally friendly

Sterile individually packed, available non-sterile\*




\*Contact Aaxis sales consultant for lead time and MOQ

PRODUCT CODE	MATERIAL	STERILE REINFORCED		SIZE	COTTON / PAPER TOWEL	AAMI LEVEL	SKU/TU (MOQ)	SKU/ CTN
S71352	50gsm 5 Layer SMMS	Yes	Yes	Medium	1 cotton	Level 3	22	22
S71353	50gsm 5 Layer SMMS	Yes	Yes	Large	1 cotton	Level 3	22	22
S71354	50gsm 5 Layer SMMS	Yes	Yes	X-Large	1 cotton	Level 3	22	22
S71355	50gsm 5 Layer SMMS	Yes	Yes	XX-Large	1 cotton	Level 3	22	22

# SureSafe Poly-Reinforced Surgical Gown

## Product Dimensions

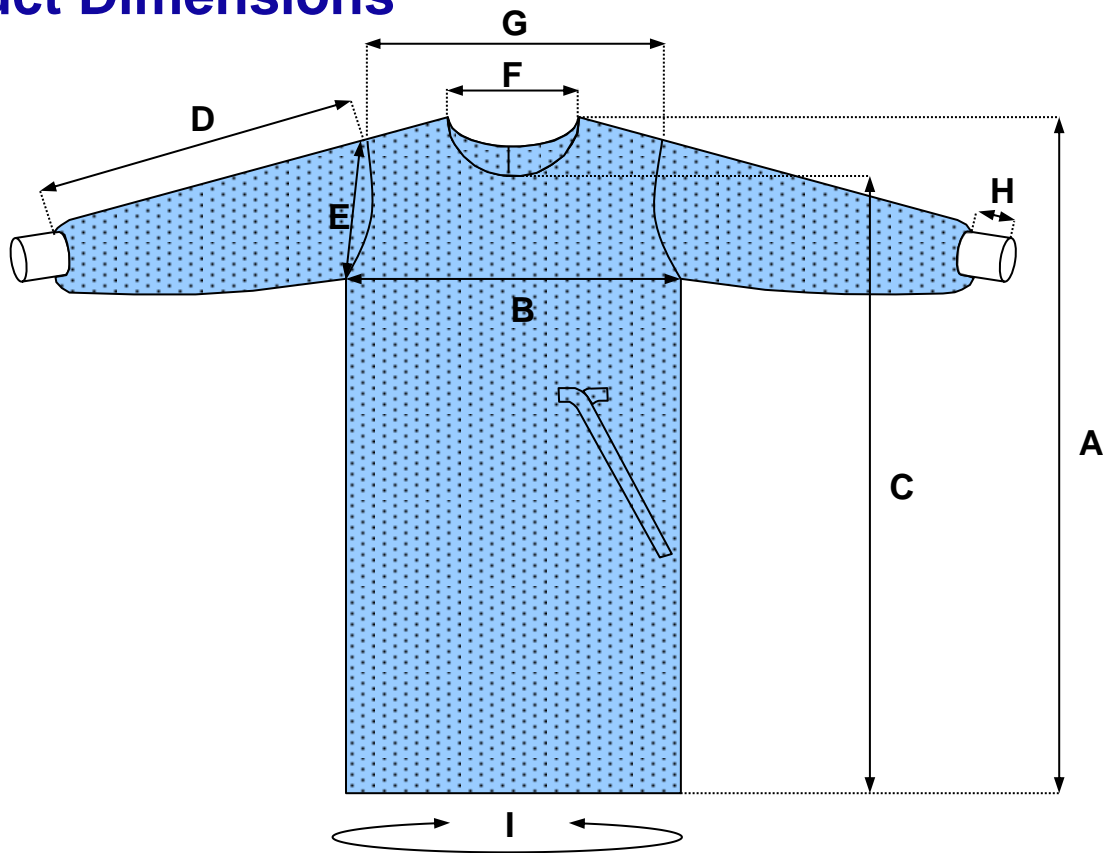

Product Code	SIZE	DIMENSIONS								
		A	B	C	D	E	F	G	H	I
S71352	Medium	120	72	112	55	26	22	57	8	150
S71353	Large	128	72	118	58	27	24	65	8	156
S71354	Extra Large	143	72	134	60	28	24	64	8	160
S71355	Extra-Extra Large	148	85	140	65	32	26	77	8	175
	Tolerance	±2.0	±1.5	±2.0	±1.5	±1.0	±1.0	±1.5	±0.75	±2.0



# Poly-Reinforced Surgical Gown

