



SureSafe Premium Gown

SureSafe® Breathable Viral Barrier Surgical Gowns offer the highest level of protection with breathable monolithic film.

BVB surgical gowns are made from BVB tri-laminate material, combining strength with softness, high breathability and excellent barrier properties against real world viruses causing infectious diseases such as HIV/AIDS (human immunodeficiency virus), SARS (Severe Acute Respiratory Syndrome) and Avian Influenza. The outer layer provides additional fluid repellence and strength, the inner layer provides a soft touch to the wearer's skin.

The layer between the inner and outer layers is a membrane which has a monolithic, non-porous structure that provides an impervious barrier, thus blocking the passage of viruses and bacteria.

The structure of the film allows moisture vapour to pass through, allowing the wearer to remain comfortable whilst providing the highest level of protection.

BVB gowns feature sonically-sealed seams (instead of stitched seams) to prevent strike through where the sleeves are joined to the body of the gown.

Tri-laminate - SPP/ polyester copolymer film / SPP

Compliance with ASTM F1671 – Maximum performance test

Compliance AAMI Level 4 - EN 13795 – Full Coverage

Abrasion and puncture resistant

Cool and comfortable

Impervious barrier

Lightweight, soft and strong

Sonically-sealed seams prevent strike through

Virtually lint free

Velcro neck binding

Book fold

Large transfer wrap – 90 x 90cm

Sterile individually packed, available non-sterile*


*Contact Aaxis sales consultant for lead time and MOQ



| PRODUCT CODE | MATERIAL | REINFORCED | SIZE | PAPER / COTTON TOWEL | AAMI LEVEL Full Coverage | SKU/TU (MOQ) | SKU/ CARTON |
|--------------|--------------------|------------|----------|----------------------|--------------------------|--------------|-------------|
| S71952 | 69gsm Tri-Laminate | N/A | Medium | 1 cotton | Level 4 | 24 | 24 |
| S71953 | 69gsm Tri-Laminate | N/A | Large | 1 cotton | Level 4 | 24 | 24 |
| S71954 | 69gsm Tri-Laminate | N/A | X-Large | 1 cotton | Level 4 | 24 | 24 |
| S71955 | 69gsm Tri-Laminate | N/A | XX-Large | 1 cotton | Level 4 | 24 | 24 |

SureSafe Premium Gown

Product Dimensions

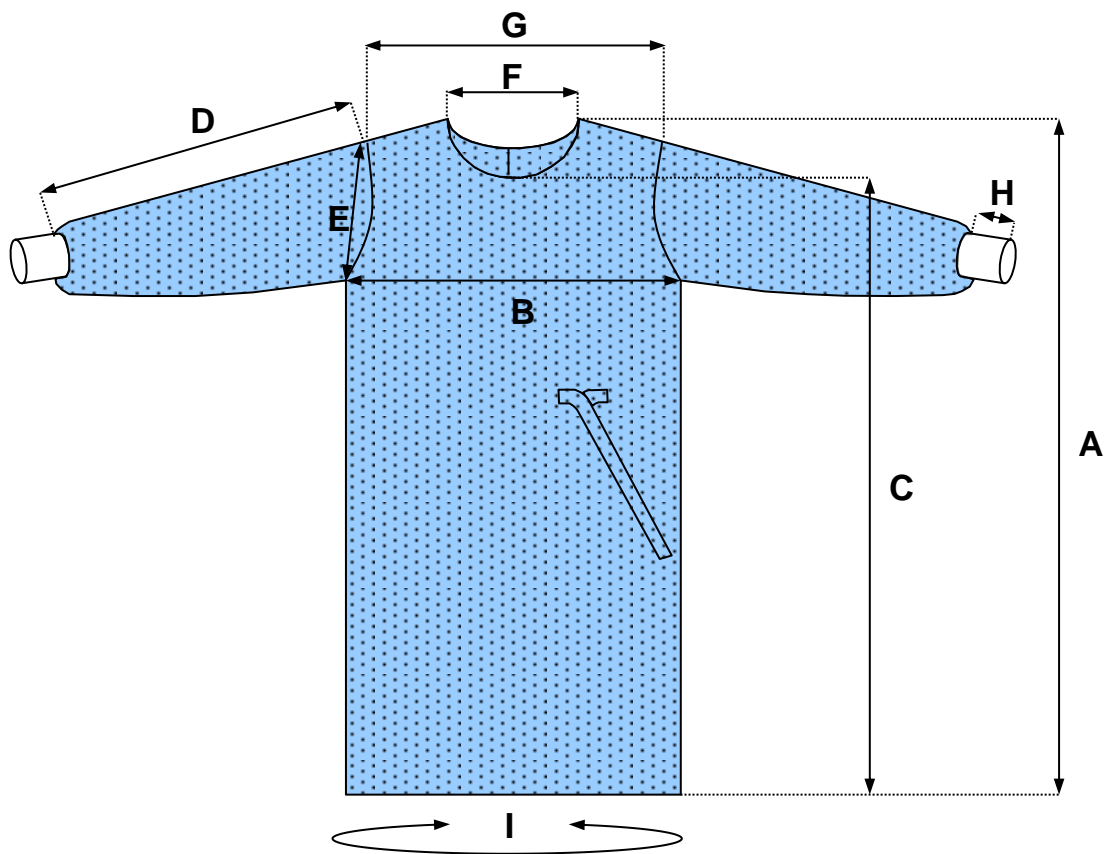

| Product Code | SIZE | DIMENSIONS | | | | | | | | |
|--------------|-------------------|------------|------|------|------|------|------|------|-------|------|
| | | A | B | C | D | E | F | G | H | I |
| S71952 | Medium | 120 | 72 | 112 | 55 | 26 | 22 | 57 | 8 | 150 |
| S71953 | Large | 128 | 72 | 118 | 58 | 27 | 24 | 65 | 8 | 156 |
| S71954 | Extra Large | 143 | 72 | 134 | 60 | 28 | 24 | 64 | 8 | 160 |
| S71955 | Extra-Extra Large | 148 | 85 | 140 | 65 | 32 | 26 | 77 | 8 | 175 |
| | Tolerance | ±2.0 | ±1.5 | ±2.0 | ±1.5 | ±1.0 | ±1.0 | ±1.5 | ±0.75 | ±2.0 |



SureSafe Premium Gown

