

Asguard®

Hydroporous Foam *(Instructions For Use)*


Asguard Hydroporous Foam dressing without border (Type I) is a highly absorbent dressing. It consists of polyurethane foam and a waterproof backing film with high permeability to allow for an effective exudate management of moderate to highly exuding wounds, and provides a moist environment to promote wound healing. It can reduce pressure on the wound during use.

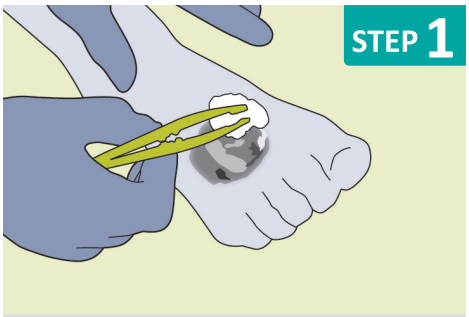


Clean the wound area, hands and instruments following normal procedures. Ensure the peri-wound skin is dry.

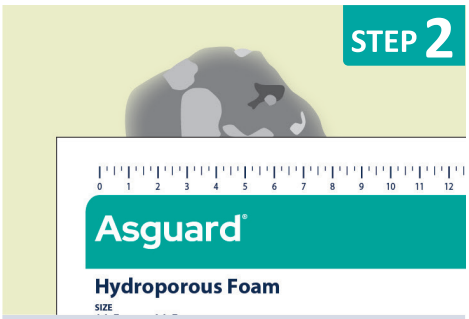


Select a size of foam dressing large enough to cover the wound edges at least 2 - 3 cm. If necessary, cut to shape.

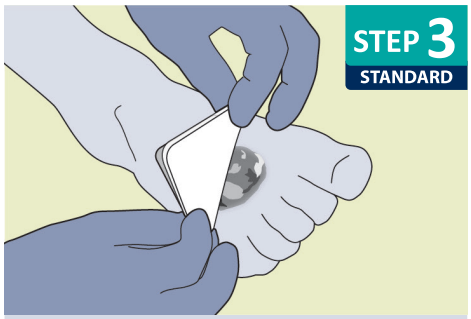


Gently apply the foam side (white side or yellowish side) to the wound ensuring good contact and fix it.

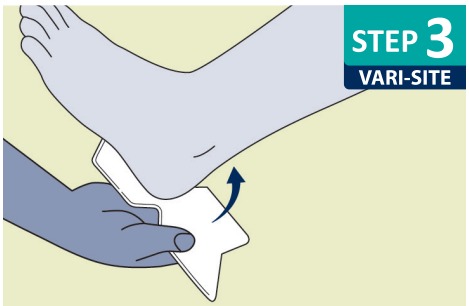

Position the narrowest part of the dressing at the base of the heel (or joint) with the white side toward the wound.

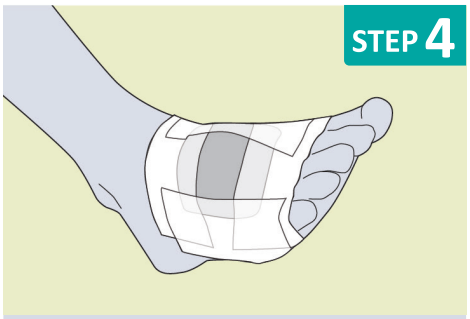

Hydroporous Foam dressing may be fixed by any secondary dressing or tape, according to clinical preference, for example: roll gauze, elastic or cohesive bandage, film dressing, tubular bandage or surgical tape.

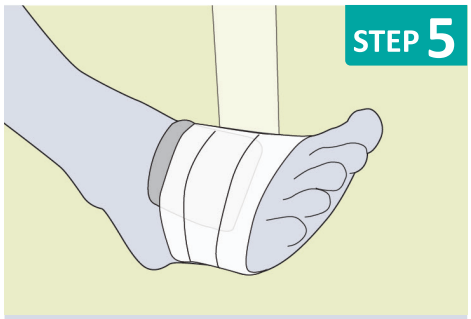

To remove the dressing: First, remove the secondary dressing, and then slowly remove the foam dressing.

CHANGE FREQUENCY

The dressing change interval may be several days. Change the dressing before it is fully saturated, at signs of leakage, or as indicated by clinical practice. A single piece dressing cannot be used for more than 7 days, and the foam dressing's accumulated use time cannot be more than 30 days.

STORAGE & TRANSPORT

The product should be stored in dry conditions.

Keep the product away from direct sunlight and keep dry.

DISPOSAL

Disposal should be handled according to local environment procedures.

At Sentry, we are dedicated to upholding sustainable principles to safeguard our planet, people, and future, and we strive to consistently act responsibly in doing so.

POUCH - PAPER	POUCH - FILM	BOX
RECYCLE ME	GENERAL WASTE	RECYCLE ME

Winner Medical Co., Ltd.
Winner Industrial Park,
No.660 Bulong Road, Longhua District,
518109 Shenzhen, China

AU REP **Sentry Medical Pty Ltd**
Suite 303, 62 Norwest Boulevard, Norwest NSW 2153, Australia
1300 995 999 | sales@sentrymedical.com.au
www.sentrymedical.com.au

Sentry MEDICAL

SMHO-MKT-IFU-0505_REV0_05_26 Seek advice from a health professional before use. All product information is correct at date of publication. The image may differ slightly from the actual product.

Asguard®

Hydroporous Foam *(Instructions For Use)*

FUNCTION PRINCIPLE

Asguard Hydroporous Foam dressing can effectively absorb wound exudate and provides a moist environment to promote wound healing. A highly breathable top film helps to reduce the risk of maceration of the wound and peri-wound as well as provides an effective barrier to bacteria to reduce the risk of infection.

INTENDED PURPOSE & INDICATIONS

INTENDED PURPOSE

Asguard Hydroporous Foam dressing is a wound dressing that absorbs wound exudate and provides a moist wound healing environment to promote wound healing.

INDICATIONS

Asguard Hydroporous Foam dressing without border (Type I) is designed for the management of moderate to heavily exuding wounds such as leg ulcers and pressure ulcers.

INTENDED USER

The product shall be used by or under the supervision of a healthcare professional.

INTENDED PATIENT POPULATION

Patient who has the indications: moderate to heavily exuding wounds such as leg ulcers and pressure ulcers.

CLINICAL BENEFIT

- Absorbs wound exudate and provides a moist wound healing environment to promote wound healing.
- Provide patient better comfort and reduce damage to the wound or surrounding skin during dressing application or removal.

CONTRAINDICATIONS

- Do not use it in the event of symptomatic fungal infection.
- Do not use it on patients with a known hypersensitivity to the product itself or its components.
- Do not use on third degree burns.

COMPLICATIONS OR SIDE-EFFECTS

Very occasionally, there may be some potential risks associated with the use of foam dressing, such as the following:

- Adverse reactions (eczema, erythema, erosion, blister, skin rash, hypersensitivity, etc.)
- Wound infections
- Risk of maceration (fluid leakage)
- Risk of skin stripping

PRECAUTIONS / WARNINGS

- Do not reuse. Reuse will cause cross-contamination.
- Do not use if inner package is damaged or opened prior to use. Do not re-sterilize.
- Do not use it together with oxidizing agents such as hypochlorite solutions or hydrogen peroxide.
- Do not use after expiry date. If the product is used after the expiry date product properties cannot be ensured.
- In case of signs of infection, consult a health care professional for adequate infection treatment.

OTHER INFORMATION

If any serious incident occurs in relation to the device, please report it to the manufacturer and the competent authority.

At Sentry, we are dedicated to upholding sustainable principles to safeguard our planet, people, and future, and we strive to consistently act responsibly in doing so.

POUCH - PAPER	POUCH - FILM	BOX
RECYCLE ME	GENERAL WASTE	RECYCLE ME

Winner Medical Co., Ltd.
Winner Industrial Park,
No.660 Bulong Road, Longhua District,
518109 Shenzhen, China

AU REP **Sentry Medical Pty Ltd**
Suite 303, 62 Norwest Boulevard, Norwest NSW 2153, Australia
1300 995 999 | sales@sentrymedical.com.au
www.sentrymedical.com.au

Sentry MEDICAL

SMHO-MKT-IFU-0493_REV0_05_26 Seek advice from a health professional before use. All product information is correct at date of publication. The image may differ slightly from the actual product.